

# DUPLEX® Drug Delivery System

Easy antibiotic therapy for compliance and patient safety.



COMING SOON:

## Cefotetan

Injection and  
Dextrose Injection in  
the DUPLEX container,  
1 gram and 2 grams

Also Available...

Cefoxitin  
Cefazolin  
CefUROXime  
CefTRIaxONE

See reverse side for  
details and ordering  
information

Ready-To-Use IV Drug Delivery System available in a portfolio covering three generations of cephalosporin products.

### Safety for Your Patient.

- Prefilled with the accurate drug and diluent doses — diluent cannot be delivered without medication
- Barcode references final admixture to help ensure correct drug
- PVC-free, DEHP-free and latex-free IV container

### Convenience and Compliance for You.

- Room temperature storage
- No mixing vials or thawing
- Ready when you are — activated at patient's bedside
- Reduces drug waste
- Joint Commission and USP Chapter <797> compliant




# B. Covered.

DUPLEX  
Drug Delivery System

**B | BRAUN**  
SHARING EXPERTISE

# DUPLEX® Drug Delivery System

## Directions for Use

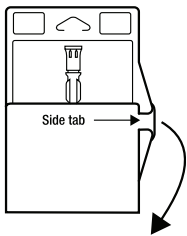

Apply patient-specific label on foil side of container. USE CARE to avoid activation. Do not cover any portion of foil strip with patient label. Unlatch side tab and unfold DUPLEX Container.

**TIP:** DUPLEX container can be stored at room temperature in Pyxis and other automated machines. If refrigerated, allow drugs to come to room temperature prior to activation.

### "PEEL"

Peel foil strip from drug chamber.

**TIP:** Due to bedside activation, if foil strip is removed and product has not been activated, product must be used within thirty (30) days for CefUROXime and CeFOXItin or seven (7) days for CefTRlaxONE and Cefazolin, but not beyond the labeled expiration date. Re-fold to protect from light and latch the side tab until ready to activate.



### "POP"

Unfold the DUPLEX container and point the set port downward. Starting at the hanger tab end, fold the DUPLEX container just below the diluent meniscus, trapping all air above the fold.

To activate, squeeze the folded diluent chamber until the seal between the diluent and powder pops open, releasing diluent into the drug powder chamber.



### "SHAKE"

Shake thoroughly to mix.

**TIP:** Visually inspect the reconstituted solution for particulate matter. Drug powder should be completely dissolved.



### "POP"

Point the set port down. Squeeze the folded DUPLEX container just below the solution meniscus, trapping all the air above the fold. This pressure will pop the second seal between the reconstituted drug and set port, releasing liquid into the port.

**TIP:** A double fold helps pop the second seal. Using aseptic technique, remove set port cover and check to assure no bag leakage prior to attaching to IV set.



### "GO"

Hang bag on IV pole.



## Ordering Information

### Cefoxitin for Injection

1 gram NDC # 0264-3123-11  
Catalog # 3123-11  
2 grams NDC # 0264-3125-11  
Catalog # 3125-11

### CefTRlaxONE for Injection and Dextrose Injection

1 gram NDC # 0264-3153-11  
Catalog # 3153-11  
2 grams NDC # 0264-3155-11  
Catalog # 3155-11

### Cefazolin for Injection USP and Dextrose Injection USP

1 gram NDC # 0264-3103-11  
Catalog # 3103-11

### CefUROXime for Injection USP and Dextrose Injection USP

1.5 grams NDC # 0264-3114-11  
Catalog # 3114-11  
750 mg NDC # 0264-3112-11  
Catalog # 3112-11

### Coming Soon!

### Cefotetan Injection and Dextrose Injection

in the DUPLEX container  
1 gram NDC # 0264-3173-11  
Catalog # 3173-11  
2 grams NDC # 0264-3175-11  
Catalog # 3175-11

All products:  
24 units per case

## Order Today!

To order or for more information call  
1-800-227-2862  
or visit  
[www.bbraunusa.com](http://www.bbraunusa.com)

**B | BRAUN**  
SHARING EXPERTISE